



# Revisiting Traditional Liquid-Liquid Extraction Techniques using SLE

Abram N. Brubaker, PhD  
MPI Research  
Mattawan, MI




# Method Development Approach at MPI Research

- Select mixed mode SPE based on functional groups present on molecule
  - For basic amine containing molecules select MCX
  - For molecules possessing a carboxyl group select MAX
  - Easier to debug a failed extraction: pH control problem most common cause
- Select SLE for the following reasons
  - Small molecular weight compound with good organic solvent solubility properties; for acidic, basic, or neutral compounds
  - SLE is preferred over HLB: partitioning differences between analyte and analog IS are minimized
  - Mixed mode or HLB extraction recoveries are extremely low due to strong interaction of analyte with the polymeric stationary phase
  - KISS principal: analysts prefer the simplicity of SLE vs SPE (automation)

## Using Ion Pairing/Back-extraction with SLE

- SLE selected because previously developed using traditional LLE
- SLE selected over SPE due to milder conditions during organic solvent removal step: time is money
- Literature example using SLE by Biotage
- Method was unsatisfactory because of analyte loss during organic solvent evaporation step
- Encountered similar problems during initial LLE

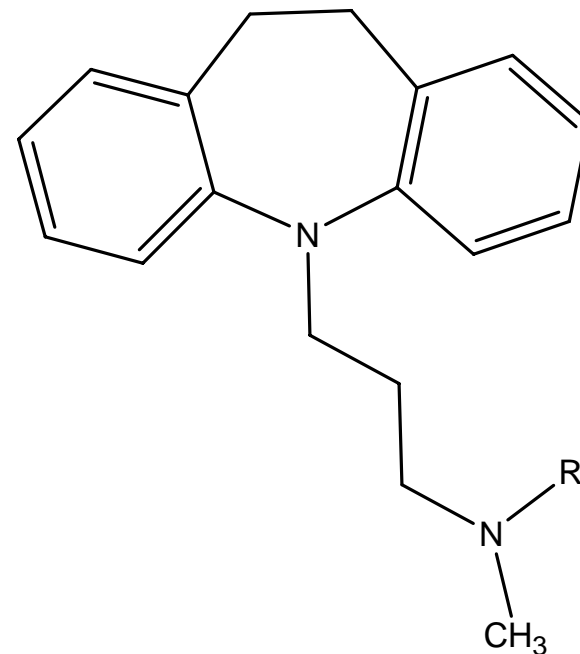


Desipramine R = H

Imipramine R = CH<sub>3</sub>



## LLE Method

- Add plasma sample to tube
- Add ion pairing agent ( $\text{ArCH}_2\text{COOH}$  in phosphate buffer pH 7)
- Add ethyl acetate
- Mix and centrifuge
- Transfer organic layer to a tube containing 5% sodium bicarbonate to back extract the ion pairing agent
- Mix and centrifuge
- Transfer organic to another tube and remove volatiles under stream of nitrogen

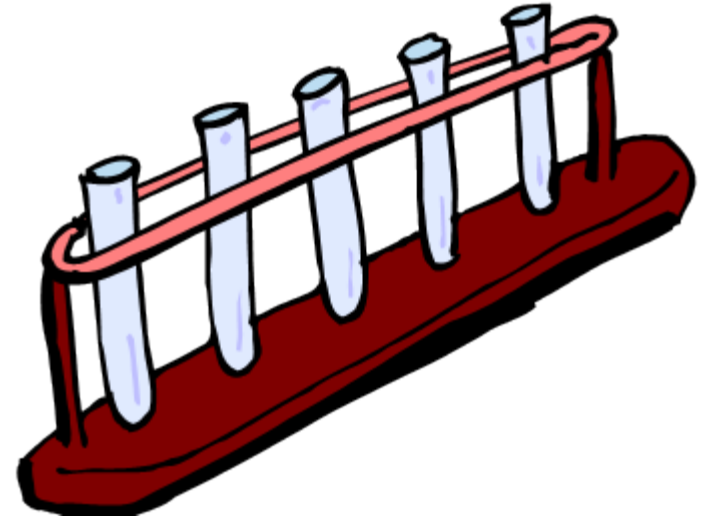

# Adapt the LLE Method to SLE

- In a 96-well plate mix
  - 50  $\mu\text{L}$  of plasma sample
  - 50  $\mu\text{L}$  internal standard solution
  - 100  $\mu\text{L}$  Ion pairing agent in phosphate buffer pH 7
- Mix and transfer to SLE plate
- Add ethyl acetate and collect organic solvent in receiver plate
- Add aqueous citric acid to receiver plate and mix well to back extract analytes
- Transfer an aliquot of lower layer to injection plate for analysis by LCMSMS






	Desipramine						Imipramine						
	QC Low 0.75 ng/mL	Rec. Std.	QC Mid 10 ng/mL	Rec. Std.	QC High 85 ng/mL	Rec. Std.		QC Low 7.5 ng/mL	Rec. Std.	QC Mid 100 ng/mL	Rec. Std.	QC High 850 ng/mL	Rec. Std.
	5536	11878	84460	173112	562849	125987 0		24766	49340	38032 3	707766	247789 0	446270 0
	4795	10867	81355	165873	616081	121629 0		24913	47254	36797 1	688349	266987 0	442743 0
	5160	11086	84862	168410	555441	121269 0		25554	46467	40075 6	691914	241937 0	439335 0
	6310	11126	82562	166537	584384	123857 0		27532	47370	38495 7	677355	252893 0	448989 0
	6050	11476	76994	161214	646029	121643 0		28456	45815	35469 7	652094	269175 0	440671 0
	6032	10802	79504	154080	625255	121778 0		27367	46508	36199 7	628213	264855 0	440255 0
Mean	5647	11206	81623	164871	598340	122693 8		26431	47126	37511 7	674282	257272 7	443043 8
SD	587	406	3014	6543	36354	18604		1551	1226	16861	29169	112945	38212
%RSD	10.4	3.6	3.7	4.0	6.1	1.5		5.9	2.6	4.5	4.3	4.4	0.9
%REC	50		50		49			56		56		58	

# SLE for Whole Blood Assay of an Unstable Compound

- Compound undergoes rapid degradation in plasma via 1,4-addition reaction
- Surprisingly stable in whole blood – stable for several weeks stored at 4° C
- Developed a LLE extraction but converted to an SLE method for automation
- Mixed whole blood, internal standard solution, and acetone
- Samples did not clog the frit even though significant residue collected on frit
- Recovery was low but greater than a protein precipitation method

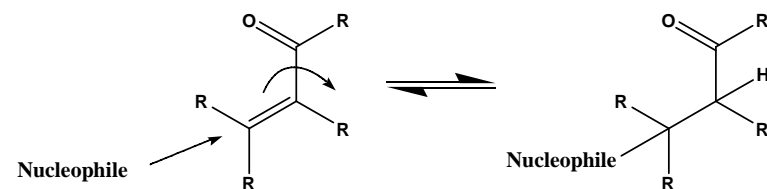

# SLE Procedure in Whole Blood


- To a 96-well plate add
  - 100  $\mu\text{L}$  of whole blood
  - 20  $\mu\text{L}$  of internal standard solution
  - 300  $\mu\text{L}$  of acetone
- Vortex, centrifuge, and transfer 200  $\mu\text{L}$  to SLE plate well and load sample onto bed
- Extract analytes by addition of 1 mL of MTBE
- Evaporate volatiles under stream of nitrogen and reconstitute sample for analysis by LC-MS/MS






	QC Low	Normalized Response	Rec. Std.		QC Mid	Normalized Response	Rec. Std.		QC High	Normalized Response	Rec. Std.
	6759	0.027	57387		151950	0.287	619462		1487230	2.93	3717320
	6196	0.024	60583		112824	0.299	525874		1666630	3.02	3555570
	9432	0.026	61213		114643	0.307	332017		1931190	3.02	3649060
	10582	0.026	57796		134994	0.302	356053		1600150	2.98	4350980
	10415	0.025	51942		123065	0.313	362181		1546860	3.00	4537660
	11688	0.027	48775		135261	0.309	523673		1696260	3.21	4872520
<b>Mean</b>	<b>9179</b>	<b>0.026</b>	<b>56283</b>		<b>128790</b>	<b>0.303</b>	<b>453210</b>		<b>1654720</b>	<b>3.02</b>	<b>4113852</b>
<b>SD</b>	<b>2219</b>	<b>0.001</b>	<b>4930</b>		<b>14863</b>	<b>0.009</b>	<b>118573</b>		<b>155543</b>	<b>0.10</b>	<b>547059</b>
<b>%RSD</b>	<b>24.2</b>	<b>4.3</b>	<b>8.8</b>		<b>11.5</b>	<b>3.0</b>	<b>26.2</b>		<b>9.4</b>	<b>3.2</b>	<b>13.3</b>
<b>%REC</b>	<b>16</b>				<b>28</b>				<b>40</b>		

# Using SLE when SPE Failed for a Complex Molecule

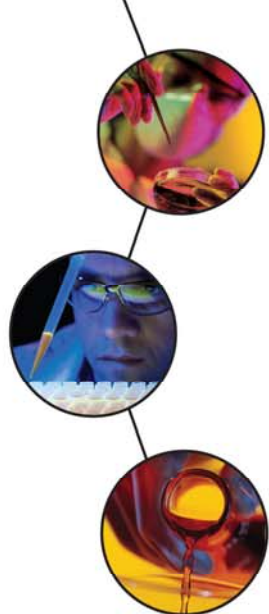

- Very polar molecule containing seven nitrogens – three basic nitrogens pKa 3 – 9
- Two aromatic rings give some lipophilic character to the molecule
- Recovery from HLB or MCX polymer based beds was very poor <50%
- Could not elute the compound even under drastic conditions – polar aprotic solvents and high concentration buffers

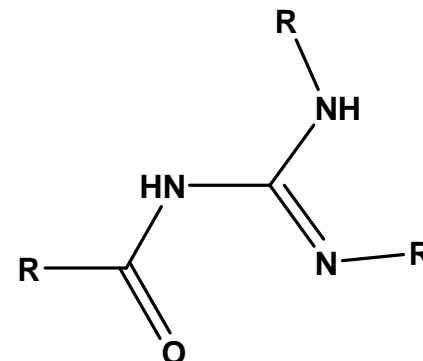




## SLE Method using Ion-Paired Extraction

- To a 96-well plate add
  - 50  $\mu\text{L}$  plasma sample
  - 50  $\mu\text{L}$  of internal standard working solution
  - 100  $\mu\text{L}$  of 0.5%  $\text{ArCH}_2\text{COOH}$  in phosphate buffer pH 7
- Mix and transfer entire volume to SLE plate
- Extract analytes with 2 x 500  $\mu\text{L}$  of MTBE
- Remove volatiles under a stream of nitrogen
- Reconstitute dried residue and analyze by LC-MS/MS





	QC Low 6 ng/mL	Low Rec. Std.	QC Mid 60 ng/mL	Mid Rec. Std.	QC High 800 ng/mL	High Rec. Std.
	26969	27579	271169	244481	3069180	3577920
	25926	27908	266247	255005	3218410	3483490
	26389	24971	268216	248272	3025840	3760070
	25887	27680	267646	257881	3068660	3744700
	28366	27573	261687	249272	3094740	3516140
	25418	28586	261662	253818	3179210	3600270
<b>Mean</b>	26493	27383	266105	251455	3109340	3613765
<b>SD</b>	1057	1241	3788	4960	73799	115342
<b>%RSD</b>	4.0	4.5	1.4	2.0	2.4	3.2
<b>%REC</b>	97		106		86	
<b>n</b>	6	6	6	6	6	6

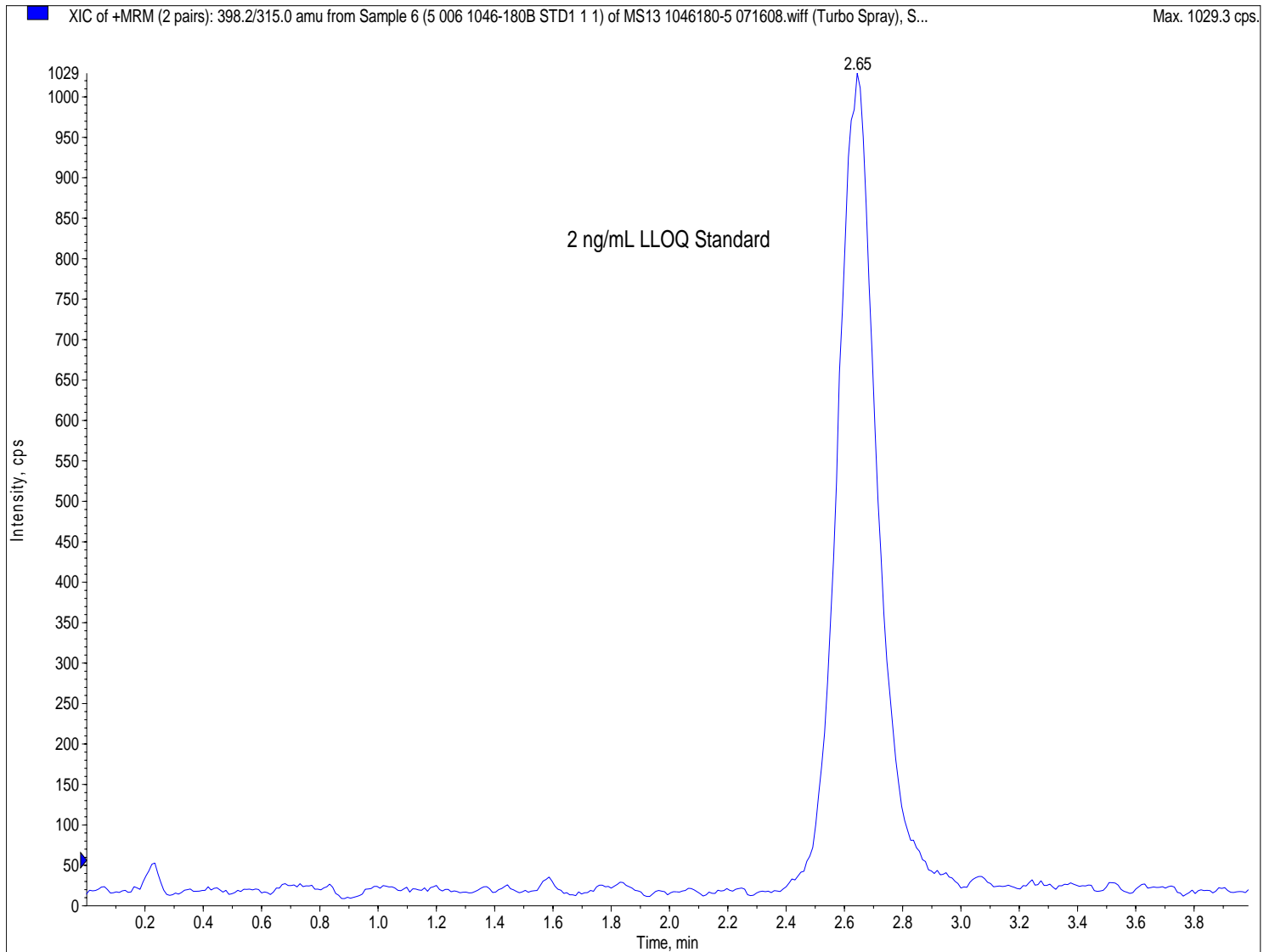

# Using SLE to Develop a Mixed Mode SPE Method

- Previous method development project using an anion exchange SPE plate failed with a phenytoin analogue
- New compound possesses active moiety similar to saccharin and should ionize at pH 4 to 5
- Chose SLE to evaluate extraction from aqueous standards to avoid secondary “chromatographic” behavioral properties of molecule and evaluate true partition coefficient properties
- Evaluated recovery of analyte at acidic, neutral, and basic pH using SLE

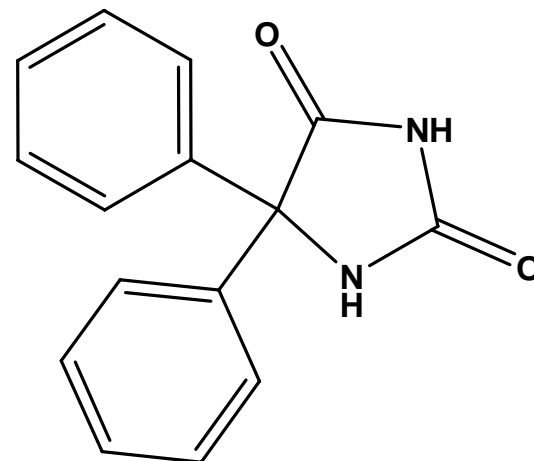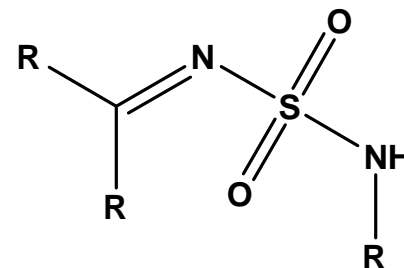



## Method:

- 1) 100  $\mu\text{L}$  of plasma (100 ng/mL analyte and IS) was diluted 1:1 with loading solvent (total volume = 200  $\mu\text{L}$ ). The entire sample was loaded onto the SLE+ plate.
- 2) Solution was allowed to flow into the bed using minimal vacuum and allowed to equilibrate for at least 5 min.
- 3) Water immiscible solvent was added to the bed and allowed to flow through. After the wells appeared dry, the solvent was allowed to equilibrate for an additional 2 min before vacuum was applied.
- 4) Samples were dried under a stream of nitrogen ( $\sim 45\text{ C}$ ) and reconstituted in 20:80 (ACN:H<sub>2</sub>O, v/v)


Acidic loading solvent = 0.1% Formic acid in Water

Neutral loading solvent = Water

Basic loading solvent = 0.5% Ammonium Hydroxide in Water



### SLE Recovery with pH Analyte

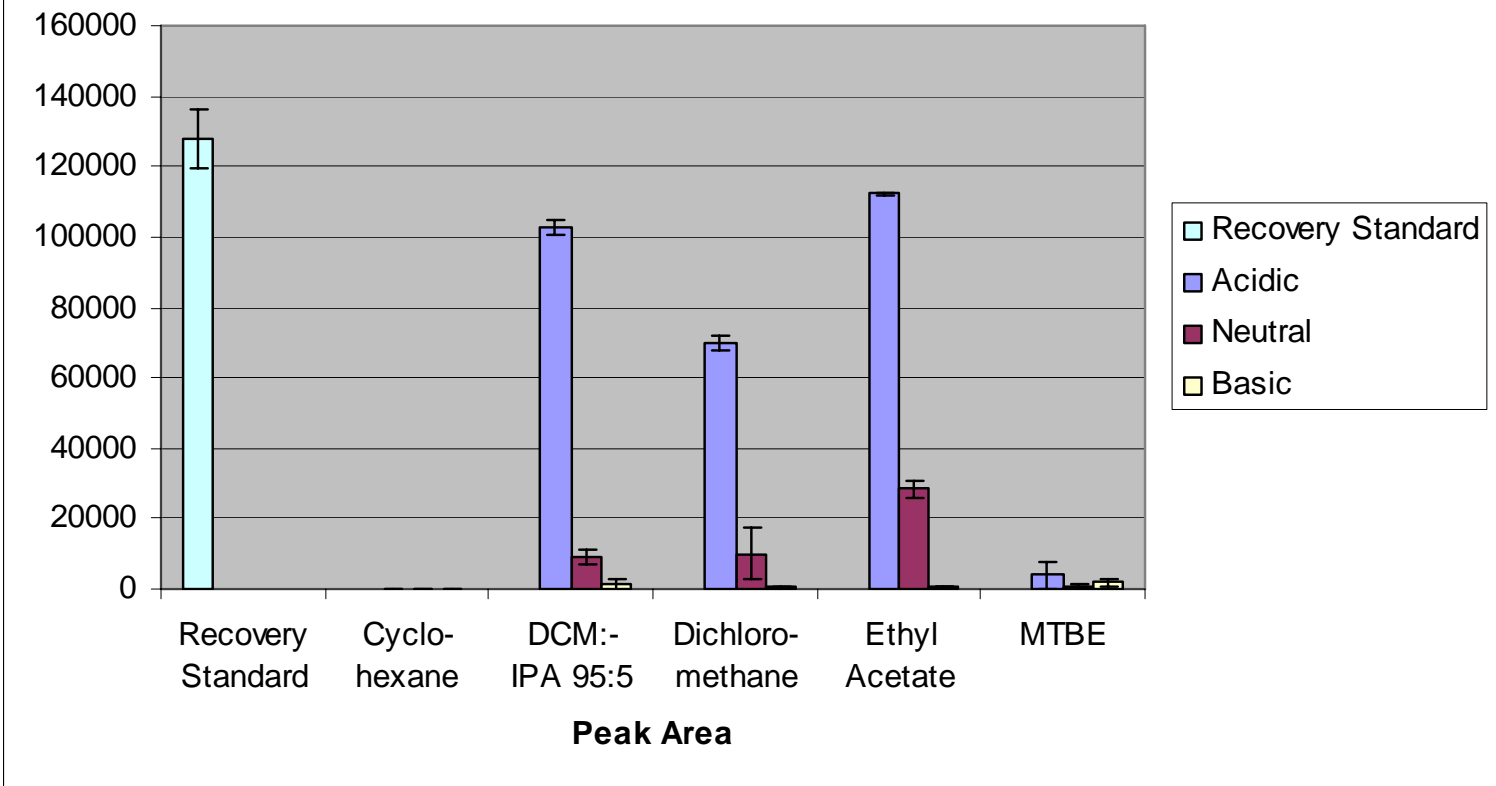





### SLE Recovery with pH Internal Standard

